

www.adventlife.net

Five-layer colourful FFP2 masks
(SERIX)

VARIETY OF COLOURS

INDIVIDUALLY PACKAGED

HUGE PRODUCTION CAPACITY



WHITE



GREEN



BLUE



BLACK

ADVENT *life*
Health in Your Hands

SERIX

Our ambition is to make healthcare widely accessible by empowering self-care globally.

Driven by our spirit of innovation, **Advent Life** specialises in self-tests for a wide range of medical conditions, while also offering other products that allow everyone to take health into their own hands. We enable people to take the first step themselves and at their own convenience.

Advent Life focuses on bringing the best from the world while prioritising safety and easy-to-use products. We design premium brands that are delivered to pharmacies, hospitals, laboratories, public organisations and government institutions.

For more information, visit www.adventlife.net.

WHAT IS COVID?

Covid-19 is an acute respiratory infectious disease caused by SARS-CoV-2. The most common symptoms are coughing and fever, but sore throat, headache and muscle aches, as well as general weakness and shortness of breath, may also occur.

Nausea, conjunctivitis, swollen lymph nodes and drowsiness, change or loss of taste and smell are some of the less common symptoms of the disease.

In some people, the infection may be asymptomatic, but it can also affect the respiratory system and other systems in the human body. In severe cases, symptoms may include shortness of breath and chest pains, and the disease may require specialist medical care or hospitalisation.

PURPOSE

Protective FFP2 masks have a design that fits tightly to the face and very effectively filters airborne particles. They are suitable for medical and daily use.

They have a traditional flared shape, with a nose clip and a double elastic ear design.

The masks in this package do not contain glass fibres, latex or other allergens. They have high bacterial filtration.

PPE INFORMATION:

Non-sterile disposable half mask made of filter material for protection against dust particles

BFE FILTER $\geq 94\%$ /dual filter

SIZE: 15.5 cm / 12 cm

COLOURS: white, blue, black, green

For medical and daily use.

Traditional convex shape, with nose clip and double elastic ear design.

Disposable, non-sterile, free of latex, glass fibres and other allergens. High bacterial filtration.

Bulgarian, English, Greek, German, Czech and Romanian languages on the packaging.

REGISTRATIONS AND CERTIFICATES

STANDARD EN 149:2001 + A1:2009

CE MARKING NOTIFIED BODY 1463

Compliance with the standard: Medical Devices Directive - 93/42/EEC EN 14683: 2019+AC:
2019 Medical face masks.

MATERIAL:



5 layers:

- 50 g of non-woven fabric
- 25 g of melted textile
- 50 g cotton with hot air
- 25 g of melted textile
- 30 g of non-woven fabric; nose clip prevents air leakage.

CONTENT:

Per unit | Hygienically packed in a blister

20 pieces in a box | Weight: 162
50 boxes in a carton | Weight: 6.5 kg
20 cartons in a pallet | Weight: 145 kg
20,000 units per pallet

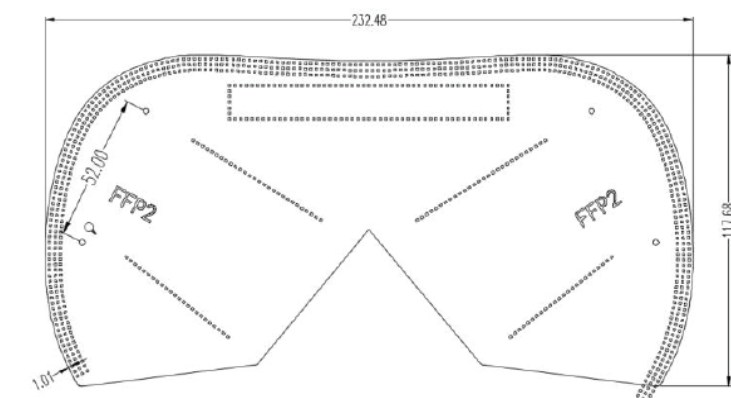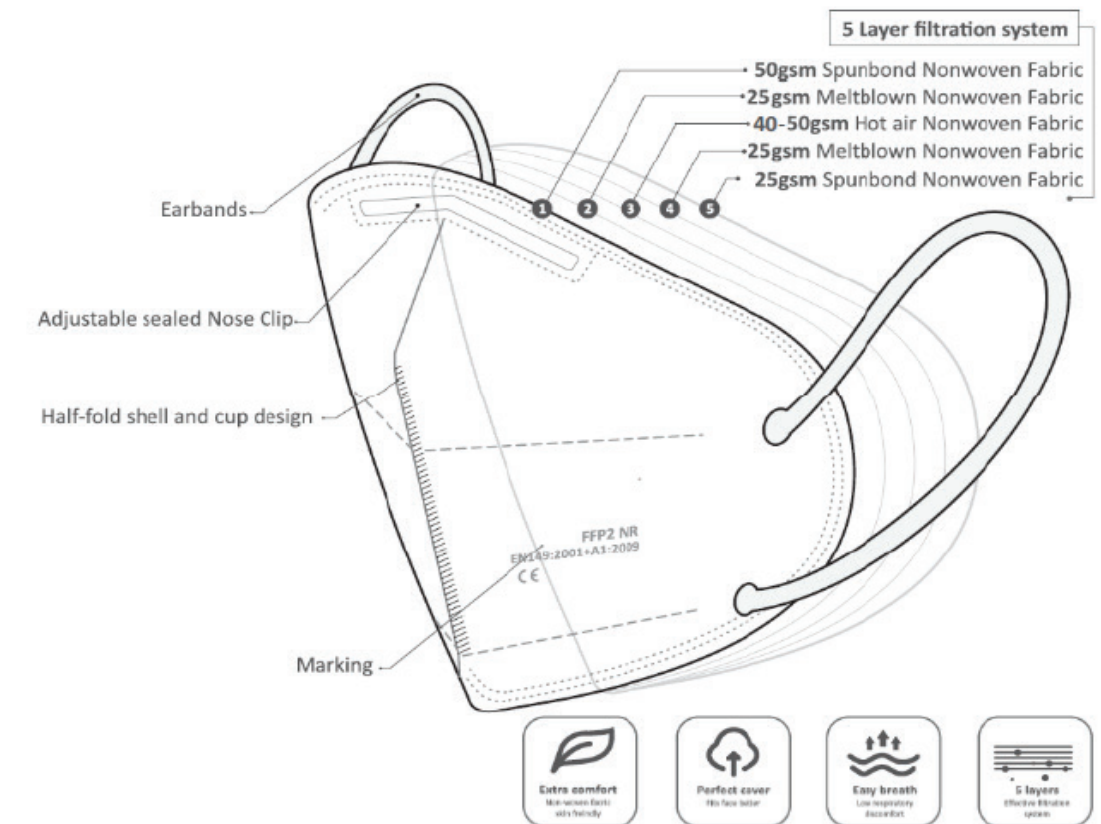

INSTRUCTIONS FOR USE:

1. Wash your hands with soap and water or a hydroalcoholic solution.
2. Open the mask and find the nose clip on top.
3. Grab and attach to the face.
4. Position it so the clip fits snugly over the nose and place the elastics behind the ears.
5. Make sure the mask is fitted correctly, fits snugly over the nose, covers the mouth and chin. After each use, dispose of the mask in a closed bin. Wash your hands with soap and water or a hydroalcoholic solution after using the mask.

STORAGE:

Store up to 25° C in a dry and ventilated place away from direct sunlight.

Cleaning/disinfection: non-sterile disposable mask. Do not clean or disinfect the mask after use.

Maximum time to wear the mask: 8 hours.

Not recommended for children under 3 years of age and people with heart and lung problems.