

VARIETY OF COLOURS

INDIVIDUALLY PACKAGED

HUGE PRODUCTION CAPACITY



GREEN



BLUE



PINK



YELLOW



BLACK

Our ambition is to make healthcare widely accessible by empowering self-care globally.

Driven by our spirit of innovation, **Advent Life** specialises in self-tests for a wide range of medical conditions, while also offering other products that allow everyone to take health into their own hands. We enable people to take the first step themselves and at their own convenience.

Advent Life focuses on bringing the best from the world while prioritising safety and easy-to-use products. We design premium brands that are delivered to pharmacies, hospitals, laboratories, public organisations and government institutions.

For more information, visit www.adventlife.net.

WHAT IS COVID?

Covid-19 is an acute respiratory infectious disease caused by SARS-CoV-2. The most common symptoms are coughing and fever, but sore throat, headache and muscle aches, as well as general weakness and shortness of breath, may also occur.

Nausea, conjunctivitis, swollen lymph nodes and drowsiness, change or loss of taste and smell are some of the less common symptoms of the disease.

In some people, the infection may be asymptomatic, but it can also affect the respiratory system and other systems in the human body. In severe cases, symptoms may include shortness of breath and chest pains, and the disease may require specialist medical care or hospitalisation.

PURPOSE

When properly fitted, a medical mask serves to block drops or splashes that may contain germs (viruses and bacteria). Thus, the mask prevents them from reaching the mouth or nose.

A medical mask also helps protect bystanders by limiting their exposure to saliva and respiratory secretions.

The medical masks in this package are hypoallergenic and do not contain glass fibres or latex. Their design includes a comfortable nose clip. The masks themselves, as well as the packaging, can be customised.

PRODUCT DESCRIPTION:

SIZE: 16.8 cm / 9.4 cm

COLOURS: blue, green, pink, yellow, black

BACTERIAL FILTRATION EFFICIENCY (BFE) (%): 99.98

INTAKE PRESSURE (KPA): 30

Nose clip

No glass fibres, no latex, hypoallergenic

English, French, German and Hungarian packaging languages

Customisation of the packaging is possible

The packaging can be printed in any language

REGISTRATIONS AND CERTIFICATES

STANDARD EN 149: 2001 + A1: 2009

In accordance with the standard: Medical Devices Directive -93/42/EEC EN 14683: 2019+AC: 2019

Medical face masks

MATERIAL:


3 layers:

- Outer: Non-woven polypropylene with spin bond
- Medium: Melt-blown non-woven polypropylene
- Inner: Non-woven polypropylene with spin bond

Adjustable nose clip.
Elastics for the ears.

CONTENTS OF 1 PACKAGE:

50 pieces in a box | Weight: 250 g
40 boxes in a carton | Weight: 10.00 kg
20 cartons in a pallet | Weight: 215 kg
40,000 pieces per pallet

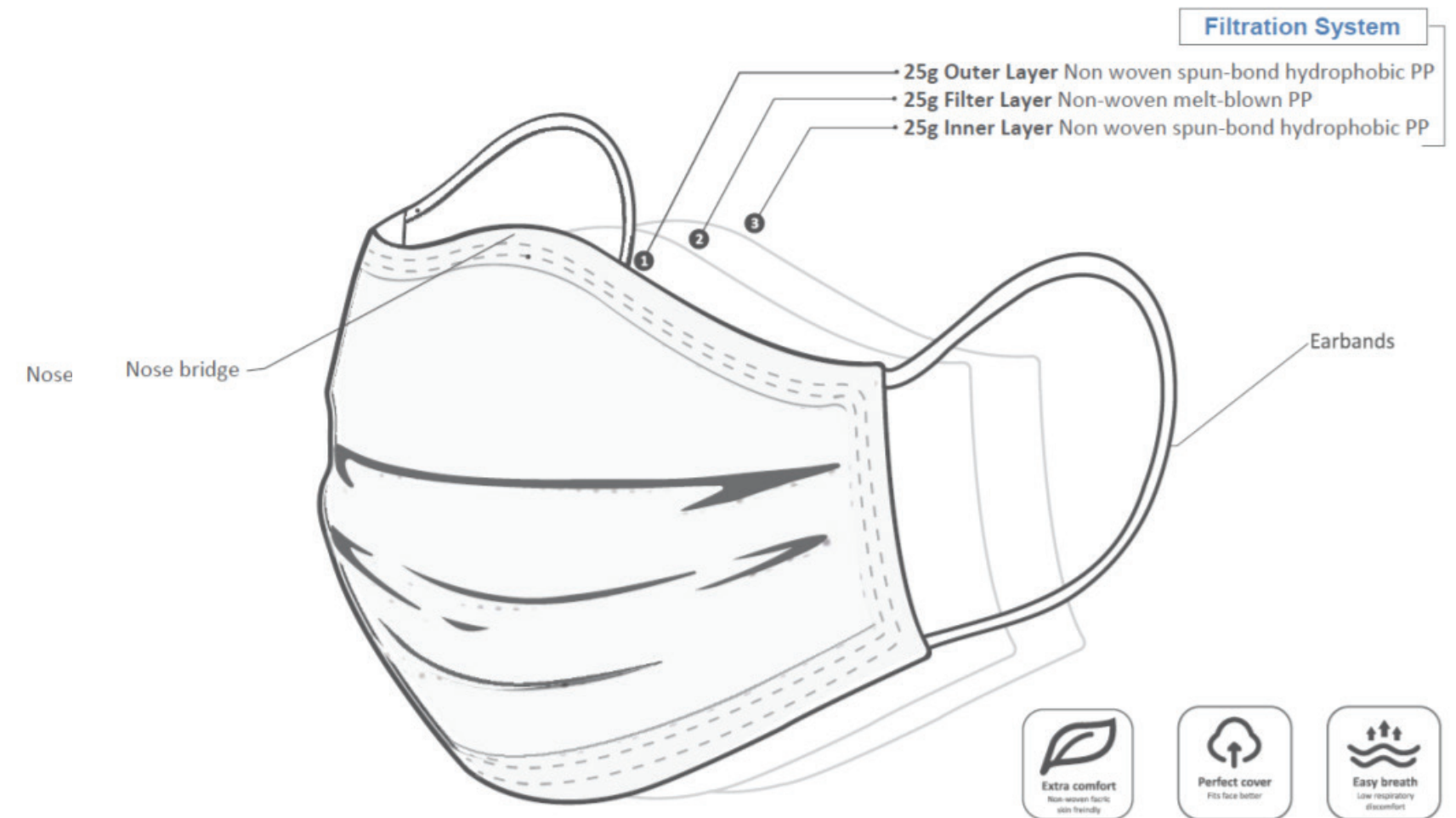

INSTRUCTIONS FOR USE:

1. Wash your hands with soap and water or a hydroalcoholic solution.
2. Adjust the elastics behind the ears.
3. Put the mask on so that it covers the nose and mouth.
4. Pull the lower part of the mask under the chin.
5. To remove the mask, grasp the elastics behind the ears (do not touch the front of the mask).
6. After each use, dispose of the mask in a closed bin. Wash your hands with soap and water or a hydroalcoholic solution after each use.

STORAGE:

Store up to 25 °C in a dry and ventilated place away from direct sunlight.

Cleaning/Disinfection:

Non-sterile disposable mask. Do not clean or disinfect the mask after use.

Maximum time to wear a mask: 4 hours

Not recommended for children under 3 years of age and people with heart and lung problems.